

Harp

Diagnostic Audiometer




Orchestral perfection
 A full line of Audiology Equipment

- PTA and speech audiometry in AC, BC and FF
- Internal memory for the speech mat. (Plus only)
- Supra-threshold tests
- Hy-Tech: can be operated through the PC



— Harp is an advanced diagnostic audiometer, capable of performing fast and accurate air, bone and speech audiometric exams, as well as several additional tests. Top-level features combined with ease of use and portability make Harp the ideal choice for a wide range of users, from hearing aid dispensers to audiologists, clinics and hospitals.

AVAILABLE TESTS


Harp is available in two configurations to suit different user needs:

- **Basic**
- **Plus**

Besides the standard audiometry exams, Harp also performs Hughson-Westlake automated audiometry, ABLB, SISI, Stenger, DLI, Tone Decay tests, and a 2 channels Master Hearing Aid. TEN Test and QuickSIN™ are available upon request only for the Plus model.

EMBEDDED SPEECH MATERIAL

Harp Plus features an integrated 4 GB flash memory for storing your speech material. Dedicated buttons allow you to access word lists with ease. The screen also displays the word currently being pronounced.

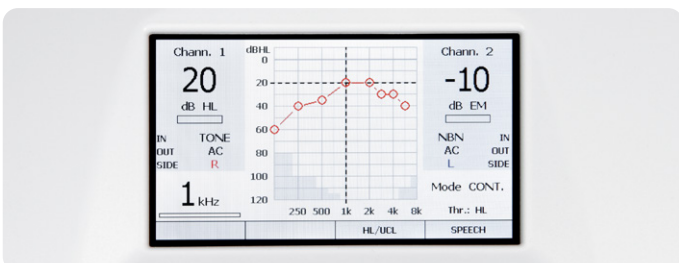

NOAH & MAESTRO INTEGRATION

Simply connect Harp to a USB port on your computer (no drivers needed), and you can transfer exam results into Maestro or Noah database. With the live view feature, you can have a wide view of the examination in progress on the computer display.

TECHNICAL DATA

AUDIOMETRY

- Intensity: AC: -10 to 120 dB HL
BC: -10 to 80 dB HL
- Freq. range: 125-8000 Hz
- Input: Tone, warble, EXT1, EXT2, Mic, flash memory (only Plus)
- Masking: NBN, WN, SN
- Output: AC L+R, BC, FF L+R, insert masking
- Talk over: Built-in or external microphone
- Monitor: Built-in speaker or headphone
- Test battery: Tone and speech audiometry, autothreshold, ABLB, SISI, Stenger, DLI, tone decay, Master Hearing Aid
only in Plus: Ten Test (opt.), QuickSIN™ (opt.)

DEVICE SPECIFICATIONS

- Data storage: Internal memory for 100 patients
- Int. flash memory: ... 4 GB (only Plus)
- Display type: Color display
- Display size/res.: ... 95 x 54 mm / 480 x 272 pixels
- Printer: Opt. integrated thermal printer
- Language: Multilingual interface
- PC interface:..... Bidirectional USB
- Modes of oper.: Stand-alone and hybrid
- Compatible sw: Inventis Maestro (stand-alone or Noah module)
- Dimensions: 32 x 32 x 9 cm
12.6 x 12.6 x 3.5 in. (DxWxH)
- Weight: without printer 1.85 kg / 4.08 lbs
with printer 2.48 kg / 5.47 lbs

Find our distributors at:
www.inventis.it/dealers

Discover the complete orchestra at: www.inventis.it

